


**Please read full instructions for use of Enerzair Breezhaler inhaler before use.** These instructions are also available by scanning the QR code or visiting: [www.breezhaler-asthma.eu/enerzair](http://www.breezhaler-asthma.eu/enerzair)

1



2

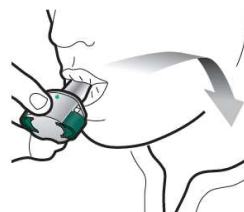

3



## Check




A hand is shown pulling the cap off the top of the device. The cap is white with a green band and a green star. The device is green and white.

You should hear a noise as the capsule is pierced. Only pierce the capsule once.



Empty

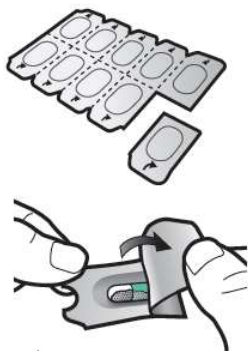

Step 1c:

**Remove capsule**

Separate one of the blisters from the blister card.

Peel open the blister and remove the capsule.

Do not push the capsule through the foil.

Do not swallow the capsule.

Breathe in quickly and as deeply as you can. During inhalation you will hear a whirring noise. You may taste the medicine as you inhale.



Step 3c:

**Hold breath**

Hold your breath for up to 5 seconds.

Step 3d:

**Rinse mouth**

Rinse your mouth with water after each dose and spit it out.

**Remove empty capsule**

Put the empty capsule in your household waste.

Close the inhaler and replace the cap.