
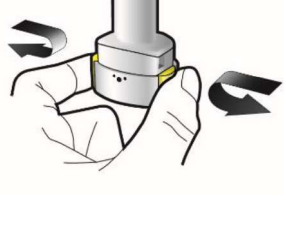
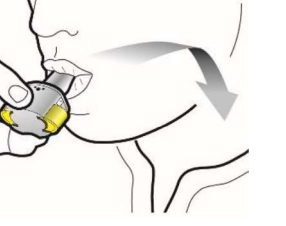
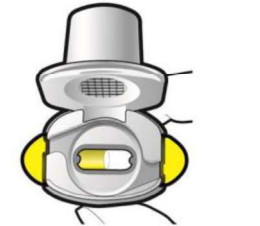
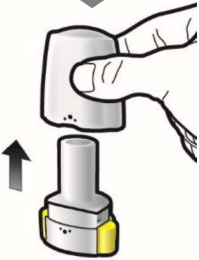
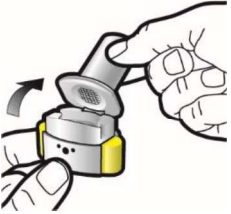
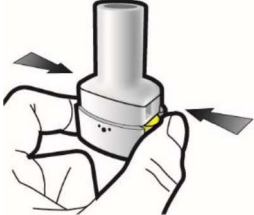



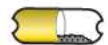
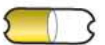


Please read the full **Instructions for Use** before using the Ultibro Breezhaler.

|  |  |  |  |
|--|--|---|---|
| Insert | Pierce and release | Inhale deeply | Check capsule is empty |
| 1 | 2 | 3 | Check |
|  <p>Step 1a: Pull off cap</p>  <p>Step 1b: Open inhaler</p> |  <p>Step 2a: Pierce capsule once Hold the inhaler upright. Pierce capsule by firmly pressing both side buttons at the same time. You should hear a noise as the capsule is pierced. <u>Only pierce the capsule once.</u></p>  <p>Step 2b: Release side buttons</p> | <p>Step 3a: Breathe out fully <u>Do not blow into the inhaler.</u></p>  <p>Step 3b: Inhale medicine deeply Hold the inhaler as shown in the picture. Place the mouthpiece in your mouth and close your lips firmly around it. <u>Do not press the side buttons.</u></p> |  <p>Check capsule is empty Open the inhaler to see if any powder is left in the capsule.</p> <p>If there is powder left in the capsule:</p> <ul style="list-style-type: none"> • Close the inhaler. • Repeat steps 3a to 3c. <div style="display: flex; justify-content: space-around;"> <div data-bbox="1134 1373 1249 1496">  <p>Powder remaining</p> </div> <div data-bbox="1321 1373 1425 1462">  <p>Empty</p> </div> </div> |



Step 1c:


Remove capsule

Separate one of the blisters from the blister card.

Peel open the blister and remove the capsule.

Do not push the capsule through the foil.

Do not swallow the capsule.

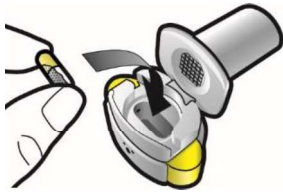

Step 1d:

Insert capsule

Never place a capsule directly into the mouthpiece.



Step 1e:

Close inhaler

Breathe in quickly and as deeply as you can. During inhalation you will hear a whirring noise. You may taste the medicine as you inhale.



Step 3c:

Hold breath

Hold your breath for up to 5 seconds.



Remove empty capsule

Put the empty capsule in your household waste. Close the inhaler and replace the cap.

Important Information

- Ultibro Breezhaler capsules must always be stored in the blister card and only removed immediately before use.
- Do not push the capsule through the foil to remove it from the blister.
- Do not swallow the capsule.
- Do not use the Ultibro Breezhaler capsules with any other inhaler.
- Do not use the Ultibro Breezhaler inhaler to take any other capsule medicine.
- Never place the capsule into your mouth or the mouthpiece of the inhaler.
- Do not press the side buttons more than once.
- Do not blow into the mouthpiece.
- Do not press the side buttons while inhaling through the mouthpiece.
- Do not handle capsules with wet hands.
- Never wash your inhaler with water.